



Knockout!



A CRISPR/Cas Gene Targeting Lab



Contents

Getting started

At a glance	P. 03
Class time requirements	P. 04
Materials needed	P. 07
Teacher prep	P. 09
Student workstation setup	P. 14
Lab cleanup	P. 14

Student guide

Background information	P. 16
Today's lab	P. 20
Student lab protocol	P. 23
Study questions	P. 26

Instructor guide

Expected results	P. 35
Troubleshooting	P. 36
Notes on lab design	P. 37
Additional student supports	P. 38
Extension activities	P. 38
Learning goals and skills developed	P. 39
Standards alignment	P. 39

At a glance

Perform cutting-edge CRISPR/Cas gene editing! This lab gives advanced students the opportunity to perform Cas9-mediated gene targeting. Students will use the CRISPR/Cas system to disrupt a gene in bacteria and observe a phenotypic change as a result.

TECHNIQUES

Micropipetting
Bacterial transformation
Gene targeting


TOPICS

Gene editing
Biotechnology

LEVEL

Advanced high school
College

Experiment overview

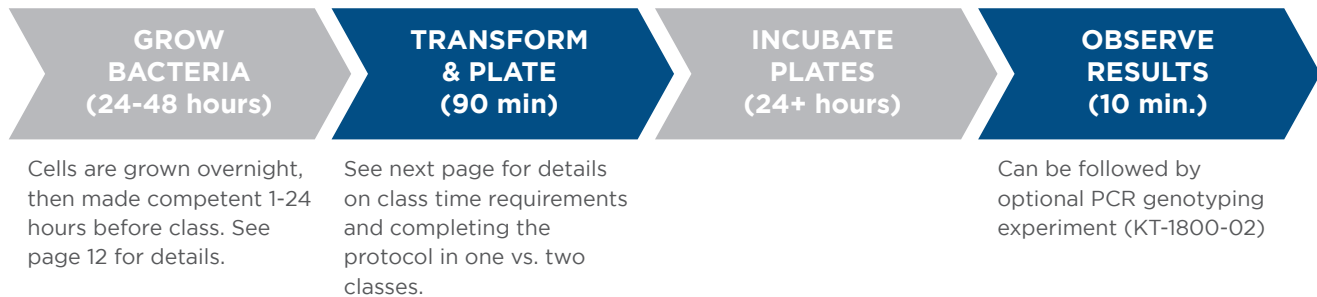


Technical support

If you have any questions about implementing this activity, contact support@minipcr.com.

Class time requirements

This protocol offers some flexibility in incubation times to help you manage the amount of class time needed for the transformation reaction. Execution of this lab is compatible with 45-minute and 90-minute class periods; please refer to the next two pages for implementation details.

Protocol steps	Time required
Add plasmid DNA to bacteria	5 minutes
Incubate on ice	30 minutes
Add recovery media	5 minutes
Recovery incubation	Minimum: 60 minutes in a standard 37 °C incubator, OR 30 minutes in a shaking 37 °C incubator Maximum: overnight at room temperature
Plate transformation reactions	10 minutes
Incubate plates at 37 °C	Minimum: 24 hours Best results after 30-48 hours
 Optional stopping point: Once colonies are visible, plates can be stored in the refrigerator for several days before students interpret the results.	
Interpret results	10 minutes


Suggested timelines for classroom implementation

FOR 90-MIN CLASS PERIODS

CLASS 1: TRANSFORMATION

Add plasmid DNA to bacteria	5 minutes
Incubate on ice	30 minutes
Add recovery media	5 minutes
Recovery incubation	30 minutes in a 37 °C shaking incubator
Plate transformation reactions	10 minutes

ALLOW COLONIES TO GROW

Incubate plates at 37 °C	Minimum: 24 hours Best results after 30-48 hours
 Optional stopping point: Once colonies are visible, plates can be stored in the refrigerator for several days before students interpret the results.	

CLASS 2: DATA COLLECTION

Interpret results	10 minutes
-------------------	------------

Continued on the next page

Suggested timelines for classroom implementation (cont.)

FOR 45-MIN CLASS PERIODS


1: SET UP TRANSFORMATION

Add plasmid DNA to bacteria	5 minutes
Incubate on ice	30 minutes
Add recovery media	5 minutes
Recovery incubation	Minimum: 1 hour at 37 °C Maximum: overnight at room temperature

2: PLATE TRANSFORMATION REACTIONS

Either after school on the same day as the transformation OR in class the day after the transformation	10 minutes
---	------------

3: ALLOW COLONIES TO GROW

Incubate plates at 37 °C	Minimum: 24 hours Best results after 30-48 hours
 Optional stopping point: Once colonies are visible, plates can be stored in the refrigerator for several days before students interpret the results.	

4: DATA COLLECTION

Interpret results	10 minutes
-------------------	------------

Materials needed

- The full lab kit (catalog no. KT-1800-01) contains all reagents and consumables for eight groups to perform the lab. A refill kit (KT-1800-03) that contains consumable reagents other than LB agar is also available. See the next page for details.
- Each kit comes with one vial of lyophilized cells. All the cells must be used on the same day after being rehydrated (see page 12). If you want to split the reagents for use with two classes that don't meet on the same day, you must order an extra vial of lyophilized bacteria (cat. no. RG-1800-02).
- The kit contains three numbered bags that require different storage conditions. Refer to the table below for details.
- Reagents must be used within six months of shipment.

	Provided	Storage	
Bag 1	ChIX Mix (chloramphenicol, IPTG, X-gal in DMSO)	1 ml	Freezer, protected from light
	pCtrl plasmid	300 µl	Freezer
	pKO plasmid	300 µl	Freezer
Bag 2	LB + ampicillin	1.7 ml	Refrigerator
	SOC recovery media	1.7 ml	Refrigerator
	Calcium chloride (CaCl ₂)	900 µl	Refrigerator
	Lyophilized DH5α <i>E. coli</i> with <i>lacZ</i> plasmid	1 vial	Refrigerator
Bag 3	LB agar powder	15 g	Room temp.
	Sterile spreaders	25	Room temp.
	Sterile Petri dishes	20	Room temp.
	Sterile transfer pipette	1	Room temp.
	Sterile 1.5 ml tubes	100	Room temp.

Required equipment

AVAILABLE AT
MINIPCR.COM

Item	Recommended quantity
Micropipettes and tips 20-200 μ l adjustable 100-1,000 μ l adjustable	1 pipette per group 1 pipette for teacher prep
37 °C incubator	Can be shared by class

Other materials supplied by user

- Distilled water
- Microwave or hot plate
- Heat-resistant flask or beaker
- Scale
- Disposable laboratory gloves
- Protective eyewear
- Heat resistant mitt
- Fine-tipped permanent marker
- 10% bleach solution

Also available

- Refill kit for CRISPR/Cas Knockout! Lab (KT-1800-03) includes cells and transformation reagents but no plasticware or solid media.
- Extra vial of lyophilized bacteria (RG-1800-02).
- Knockout! PCR Genotyping Experiment (KT-1800-02).

Teacher prep

- Reagents are sufficient for eight lab groups.
- The teacher prep for this activity takes place over several days and must start 1-2 days before the lab.
- The table below provides an overview of the teacher prep and the subsequent pages provide detailed instructions.

Prep	Time required	Timeline
Prepare plates	30 minutes	Can be completed up to one week in advance
Prepare bacteria		
A. Rehydrate bacteria	20 minutes	Must be completed the day before you plan to add CaCl ₂ to the cells
B. Add CaCl ₂	10 minutes	Minimum: 1 hour before the lab Maximum: 24 hours before the lab Longer incubation in CaCl ₂ increases transformation efficiency
Dispense reagents	10 minutes	Must be completed before class on the day of the lab


Continued on the next page



Protective gloves and eyewear should be worn for the entirety of this experiment.

Up to one week before the lab: prepare plates

- You will need two plates per lab group.
- IMPORTANT NOTE: The instructions that follow are for preparing eight LB agar plates at a time to allow the use of a smaller bottle that will fit in a standard microwave.
- TIP: If a 500 ml bottle is too tall to fit in your microwave, you can prop it at an angle by placing it inside another container.

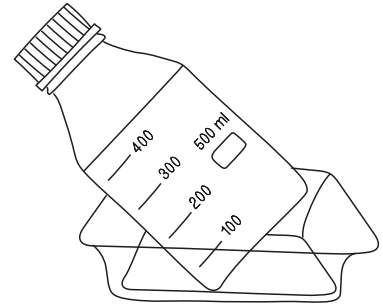

Materials needed for this section:

From Bag 1 of lab kit (stored in freezer):

- ChIX mix

From Bag 3 of lab kit (stored at room temp):

- LB agar powder
- Sterile Petri dishes

Equipment and materials supplied by user:

- Distilled water
- Microwave or hot plate
- Heat resistant mitt
- 100-1,000 μ l micropipette and tips

1. Thaw the ChIX mix, and then set the tube aside. Note that this can take some time. We recommend placing the tube in an incubator or heat block set at 37 °C.
2. Add 190 ml distilled water to a heat-proof bottle with at least a 500 ml capacity.
3. Add 7.5 g LB agar powder and swirl to mix.
4. Loosely cap the bottle to allow for steam to vent.
5. Microwave in 30-second increments just until the solution boils. Continuously monitor to ensure the solution doesn't boil over.



Caution: The solution will be very hot!

6. Using a heat resistant mitt, carefully swirl the solution.
7. Repeat steps 5 and 6 at least two more times. Continue until the LB agar powder is fully dissolved and the solution is transparent.
8. Let the LB agar solution cool until you can touch the bottle with your bare hands, but not so long that the agar begins to set. This should take approximately 3-5 minutes.

9. Mix the ChIX, then use a micropipette to add 475 μl to the LB agar solution and swirl to mix. The ChIX may refreeze at room temperature. Make sure the solution is fully thawed before using.
10. Pour LB agar into sterile Petri dishes. The LB agar in each dish should be ~3 mm thick. You should be able to pour at least 8 plates with this volume of LB agar.
11. Immediately put the lids on the Petri dishes and allow LB agar to solidify. If possible, let the plates sit at room temperature for several hours. This helps prevent excessive condensation from forming while stored in the refrigerator.
12. For a class with eight lab groups, you will need 16 plates total. Repeat steps 2-11 to pour another set of eight plates.
13. Once agar is solid, stack plates upside down with the agar on top.
14. Place stacked plates in a plastic bag and store in the refrigerator until use. Because X-gal present in the ChIX is light sensitive, protect plates from light if possible.

One or two days before class: Prepare bacteria

A. Rehydrate bacteria

Materials needed for this section

From Bag 2 of lab kit (stored in refrigerator):

- Lyophilized *E. coli*
- LB + ampicillin

Equipment and materials supplied by user:

- 20-200 μ l micropipette and tips
- 100-1,000 μ l micropipette and tips
- 37 °C incubator

From Bag 3 of lab kit (stored at room temp):

- Sterile 1.5 ml tubes

1. Uncap the vial of lyophilized bacteria by peeling away the metal seal and then pulling out the rubber stopper.
2. Add 900 μ l of LB + ampicillin liquid media to the vial.
Note: Make sure you use LB + ampicillin, not SOC recovery media (the containers look very similar).
3. Reinsert the rubber stopper to cap the vial and invert several times to mix.
4. Aliquot 50 μ l of rehydrated bacteria into 16 sterile 1.5 ml microtubes (two tubes per lab group).
To prevent the bacteria from settling at the bottom of the vial, mix the solution regularly with a micropipette while aliquoting.
5. **Incubate at 37 °C for 16-24 hours.** Do not disturb or agitate samples during the incubation period.

B. Add CaCl₂ to the bacteria

Materials needed for this section

Rehydrated bacteria from the previous step

Equipment and materials supplied by user:

From Bag 2 of lab kit (stored in refrigerator):

- CaCl₂

- 20-200 μ l micropipette and tips
- Refrigerator

1. Add 40 μ l of CaCl₂ to each of the 16 tubes of bacteria.
2. **Place tubes in the refrigerator for a minimum of 1 hour and up to 24 hours.** Longer incubation times result in increased bacterial competency and better transformation efficiency.

The day of the lab: Dispense reagents

Materials needed for this section

From Bag 1 of lab kit (stored in freezer):

- pCtrl plasmid
- pKO plasmid

From Bag 2 of lab kit (stored in refrigerator):

- SOC recovery media


From Bag 3 of lab kit (stored at room temp):

- Sterile 1.5 ml tubes

Equipment and materials supplied by user:

- 20-200 μ l micropipette and tips

1. Thaw tubes containing the DNA samples by placing them at room temperature.
2. For each lab group, dispense the DNA into labeled sterile 1.5 ml microtubes:
 - pCtrl DNA, 30 μ l
 - pKO DNA, 30 μ l
3. For each lab group, use sterile micropipette tips to dispense 120 μ l SOC recovery media into labeled sterile 1.5 ml microtubes.
4. You can store the tubes at room temperature for up to five hours before use.

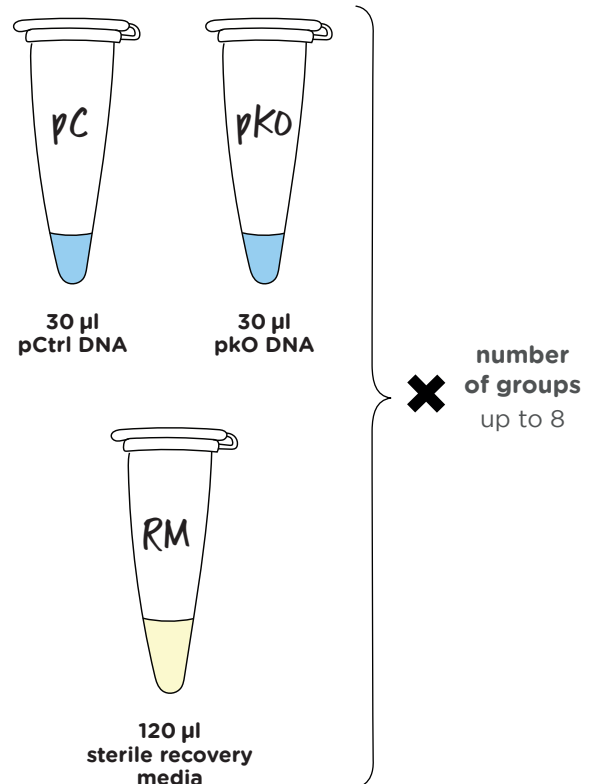

Student workstation setup

- Depending on the length of your class periods, you may choose to perform the transformation and plating steps on consecutive days.
- The lists below have been separated to make it clear what is needed for each.

To perform the transformation, every group will need:

Prepared bacteria <ul style="list-style-type: none"> • Rehydrated and grown for 16-24 hours at 37 °C • CaCl₂ added and incubated in the refrigerator for 1-24 hours 	2 tubes each with 90 µl
Plasmid DNA <ul style="list-style-type: none"> • pCtrl • pKO 	30 µl of each plasmid
SOC recovery media	120 µl in a sterile tube
20-200 µl micropipette and tips	
Fine-tipped permanent marker	
Crushed ice	
Access to a 37 °C incubator	

To plate the transformation reactions, every group will need:

ChIX LB agar plates, prewarmed	2
Sterile spreaders	2
20-200 µl micropipette and tips	
Fine-tipped permanent marker	
Access to a 37 °C incubator	

Lab cleanup

After the lab, treat bacterial plates and any materials that came into contact with bacteria in 10% bleach for at least 20 minutes.



Student guide



Background information	P. 16
Today's lab	P. 20
Student lab protocol	P. 23
Study questions pre-lab	P. 26
CER table	P. 29
Post-lab study questions	P. 31



Background information

Overview

In today's lab you will use one of biotechnology's most exciting tools: the CRISPR/Cas system. Scientists use this system to alter the DNA in living cells. Become a genetic engineer yourself when you use CRISPR/Cas to disable a gene in bacteria!

Genome editing

Deliberately altering an organism's DNA has been a major goal of scientists for decades. Scientists want to edit DNA for two main reasons. First, to understand the biological processes in any organism, we need to understand how its genes work. One of the best ways to determine what a specific gene does is to observe what happens when a gene is disabled or "knocked out". Second, once we understand the function of a gene, scientists can deliberately modify the gene to change the organism in specific ways. For example, scientists could correct harmful mutations or introduce advantageous traits like making crops resistant to drought.

Genome editing is making a specific change to an organism's DNA. In the 1990's and 2000's, scientists developed a handful of tools for genome editing, but these tools had significant drawbacks. None were able meet the goal of efficiently introducing predictable and precise DNA edits to virtually any organism. Between 2012-2013, work from several scientists demonstrated a powerful new tool for genome editing that could meet that elusive goal: the CRISPR/Cas system.

Using CRISPR/Cas as a genome editing tool

Scientists did not invent the CRISPR Cas system from scratch; like most biotechnology tools, it has a natural origin. CRISPR/Cas is found in bacteria and archaea where it plays a role in immune function, providing protection from harmful viruses. Scientists have repurposed this bacterial immune system as a genome editing tool because it allows them to target specific DNA sequences with relative ease.

The CRISPR/Cas system involves two main components: a Cas nuclease and a guide RNA (Figure 1). Nucleases are enzymes that cut nucleic acids, such as DNA, like a pair of molecular scissors. There are several Cas nucleases, derived from different species of bacteria, that work in slightly different ways, but Cas9 is most commonly used in genome editing. Cas9 is a powerful tool because it can be programmed to specifically cut nearly any DNA sequence.



What's in a name?

CRISPR stands for **clustered regularly interspaced short palindromic repeats**—what a mouthful! It refers to an area of the bacterial genome involved in the immune defense against viruses. This bacterial defense mechanism relies on two main components: the DNA region we call CRISPR and the Cas9 nuclease. The CRISPR/Cas9 genome editing technique that scientists use relies heavily on Cas9, yet the CRISPR region doesn't play a role in this type of genome editing. However, saying 'CRISPR' sure is catchier than calling it Cas9 genome editing, which is probably why this term has stuck around.

The site where Cas9 will cut is determined by a short RNA molecule called a *guide RNA* (Figure 1). The guide RNA will bind to Cas9 and form a complex that scans the genome. At one end of the guide RNA is a stretch of approximately 20 bases that determines the DNA sequence that Cas9 will cut. When this region of the guide RNA sequence encounters and binds to a complementary DNA sequence, Cas9 will cut both strands of the DNA (Figure 1).

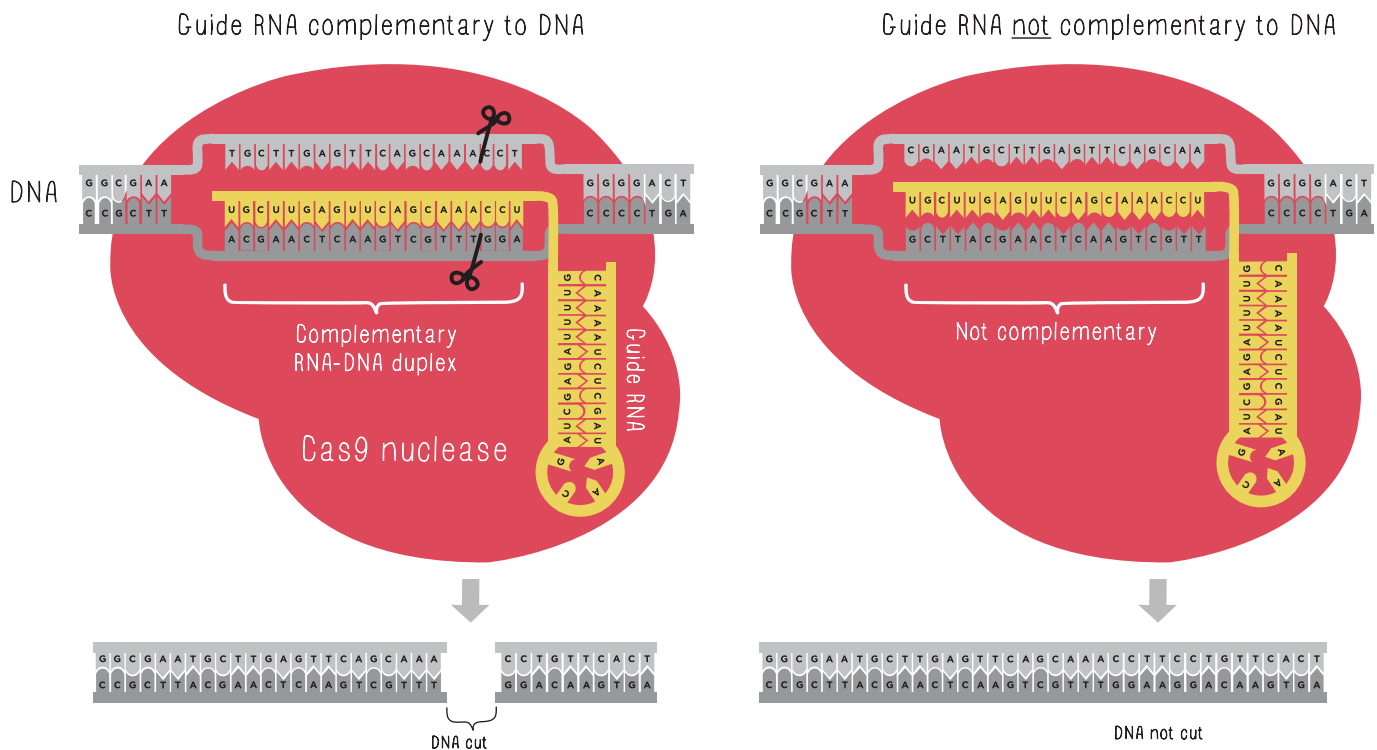


Figure 1. Cas9 cuts DNA as instructed by a guide RNA. When the guide RNA encounters a complementary DNA sequence, Cas9 will cut the DNA (left). When the guide RNA is not complementary to the DNA, Cas9 does not cut (right).

What makes CRISPR/Cas9 such a powerful genome editing tool is that scientists can customize this ~20 base guide sequence to target virtually any DNA sequence that they are interested in. Cas9 will then target that specific region of DNA, cutting exactly where the scientist planned, and not elsewhere in the genome.



To understand why cutting DNA at a specific place is important in genome editing, we have to discuss what happens in the cell after the DNA is cut. In eukaryotic cells, once Cas9 cuts the target DNA, the cell will try to repair the break (Figure 2). One way the cell can accomplish this is to reattach the broken strands of DNA through a process called non-homologous end joining (NHEJ). When the cell does this, it often ends up adding or removing a few DNA bases. These act much like typos. These typos, or mutations, can disrupt a gene or other DNA sequence.

While we often think of mutations like these as being harmful to an organism, they can sometimes be used as a tool by scientists. For example, mutations that disable a gene can help scientists understand the function of that gene. Before CRISPR/Cas, accomplishing this was a difficult and time-consuming process that could only be done in certain organisms. Because CRISPR/Cas makes it feasible to disable genes in a wide range of organisms, this is the most common use of the system.

The CRISPR/Cas system can also be used to introduce specific changes to the genome using a different type of DNA repair called homologous recombination. For more information on this topic, refer to the DNAdots article on CRISPR/Cas9 (<https://dnadots.minipcr.com/dnadots/crisprcas-9>).


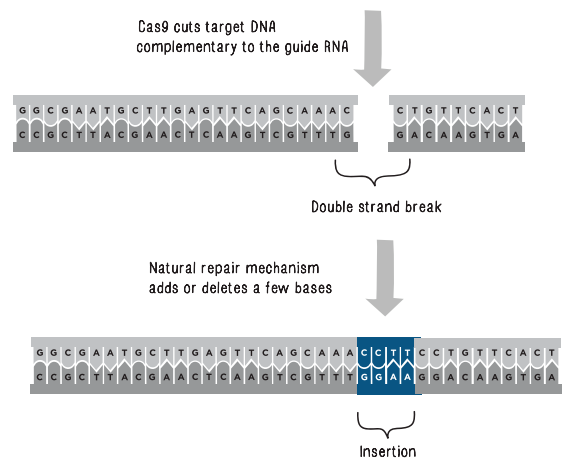


Figure 2. DNA breaks can lead to mutations.

Breaks in the DNA can be repaired by two different mechanisms in eukaryotic cells. The more common repair mechanism, non-homologous end joining, reattaches the broken strands of DNA, but in the process often adds or removes a few random DNA bases. These insertions and deletions often disable the gene.

Lefty snails: using CRISPR/Cas9 genome editing to understand gene function

Scientists had long observed that the shells of most snails spiral to the right, but a few rare snails have shells that spiral to the left. This was the case for Jeremy the garden snail, whose plight to find another lefty snail to mate with made international headlines. While scientists suspected that a gene called *Lsdia1* controlled the direction in which a snail's shell spirals, prior to CRISPR/Cas genome editing, there were no tools to edit DNA in snails. Without a way to disrupt the *Lsdia1* gene and observe the effects on shell spiraling, there was no way to test this hypothesis. Recently, scientists used CRISPR/Cas9 genome editing to knock out the *Lsdia1* gene in snail embryos. After the snails matured, they observed the direction in which their shells coiled, and definitively showed that *Lsdia1* dictates shell spiraling (Abe and Kuroda, 2019). For more information on this work, refer to the *Snail shell coiling paper model* activity at <https://www.minipcr.com/crispr-paper-model/>.

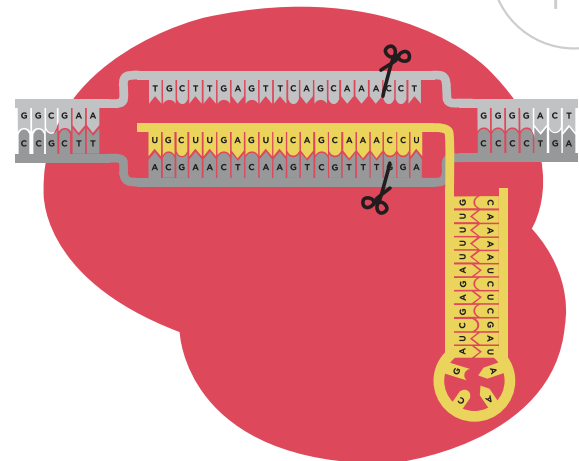

Advantages of the CRISPR/Cas system for genome editing

There are numerous ways in which the CRISPR/Cas system has revolutionized genome editing. Here we will focus on two: adaptability and specificity.

Adaptability

In theory, the CRISPR/Cas system can target any DNA sequence in any organism. Previous genome editing techniques were so complicated and inefficient that it was only practical to modify the genomes of certain organisms—this was limited to the most commonly used lab organisms such as fruit flies or mice. The application of the CRISPR/Cas system as a genome editing tool provides scientists with a relatively easy way to alter the genome of virtually any organism.

Specificity

It is essential to control with precision where genome editing tools act, to avoid introducing unwanted mutations, and until recently, this remained a substantial challenge. Because the guide RNA has a ~20 nucleotide recognition sequence, Cas9 targeting is quite specific. The chance of any specific 20 base sequence matching a random 20 base stretch is less than one in a trillion. This means that even in a genome that is billions of base pairs long, it is likely that the only place the Cas9 enzyme will cut is the specific place in the genome for which the guide RNA is designed. While the CRISPR/Cas system has high specificity compared to older genome editing technologies, it is not perfect. Sometimes a partial match between the guide RNA and the genome can lead to unexpected off-target cutting. To reduce the chance of off-target editing, scientists are experimenting with modifications to the CRISPR/Cas system to further increase its specificity.



Today's lab

Today, you will use the CRISPR/Cas9 system to disable, or knock out, a gene in *E. coli* bacteria. The *E. coli* genome contains a single circular chromosome, but this is not the only place where bacterial genes are found. Bacteria also contain plasmids, small rings of DNA that typically include at least one gene that confers an advantageous trait, such as antibiotic resistance. Bacteria can transfer plasmids to each other, sharing these beneficial traits. Scientists use plasmids as tools to introduce DNA to bacteria through a process called transformation. For more information on transformation, refer to <https://www.minipcr.com/tutorials/>.

Overview

The *E. coli* you will use already contain a plasmid called pLacZ that carries the *lacZ* gene, and this will be the gene targeted by the CRISPR/Cas system. The *lacZ* gene is a convenient target because there is a quick and easy way to tell if the gene has been disrupted. The *lacZ* gene encodes an enzyme called β -galactosidase, which catalyzes the breakdown of the sugar lactose. But β -galactosidase can also catalyze the breakdown of a chemical called X-gal. When X-gal breaks down, one of the products is blue (Figure 3). This means that the presence of blue color indicates that the *lacZ* gene is functional and β -galactosidase protein is present.

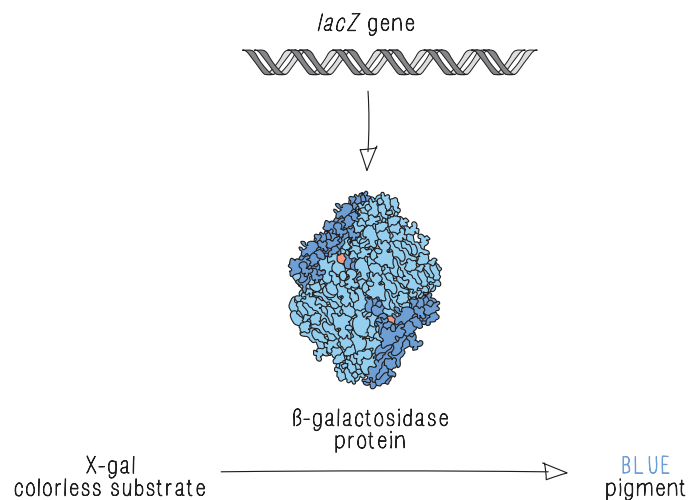


Figure 3. The *lacZ* gene. The *lacZ* gene encodes the β -galactosidase protein. β -galactosidase can process a chemical called X-gal to produce a blue pigment.

Overview of experiment

While the CRISPR/Cas system is native to bacteria, different types of bacteria use different Cas proteins. The *E. coli* bacteria you will be using today do not naturally contain Cas9. Instead, you will transform a second plasmid called pKO into the *E. coli* to introduce the *cas9* gene along with the sequence for a custom guide RNA that is complementary to 20 bases of the *lacZ* gene (Figure 4, middle). The bacteria will transcribe and translate the *cas9* gene to make Cas9 protein, while the guide sequence will just be transcribed into RNA. Then, the guide RNA will bind to Cas9 and together they will scan the DNA of the cell for a sequence complementary to the guide RNA. When the guide RNA binds to the *lacZ* sequence in the plasmid, Cas9 will cut the DNA (Figure 5, left).




You will also perform a control transformation. In this reaction, you will also transform bacteria with a plasmid called pCtrl that contains *cas9* gene and the instructions for a guide RNA with a random 20 base sequence (Figure 4, bottom). This random sequence is not expected to be complementary to any DNA in the cell, so Cas9 will not cut (Figure 5, right).

After growing the bacteria on plates containing X-gal, you will use the color of the bacterial colonies to verify whether the Cas9 protein cut the *lacZ* gene (Figure 5). In this experiment, blue colonies indicate the presence of a functional copy of the *lacZ* gene (Figure 5, right), and white colonies indicate that the *lacZ* gene has been disabled, or knocked out, by Cas9 (Figure 5, left).

Become a genetic engineer yourself. Let's knock out some genes!

Plasmid already present

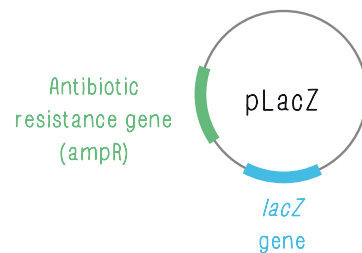

Plasmids used for transformation

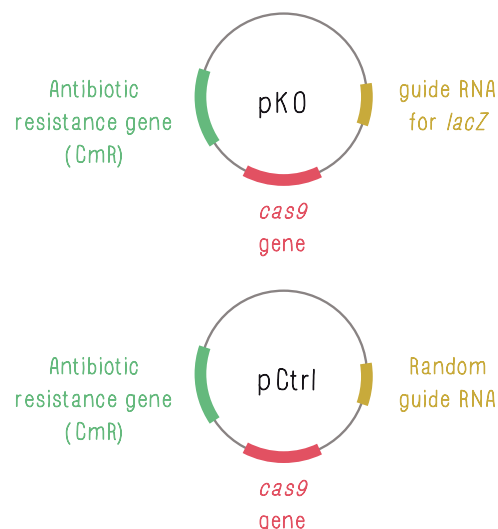


Figure 4. Plasmids used in experiment. Bacteria used in the experiment already contain the pLacZ plasmid (top) that carries the *lacZ* gene and confers resistance to the antibiotic ampicillin. These cells are then transformed with either the pKO plasmid (middle), that carries the *cas9* gene and instructions for a guide RNA complementary to the *lacZ* gene, or the pCtrl plasmid (bottom), that carries the *cas9* gene and instructions for a random guide RNA. Both the pKO and pCtrl plasmids confer resistance to the antibiotic chloramphenicol.

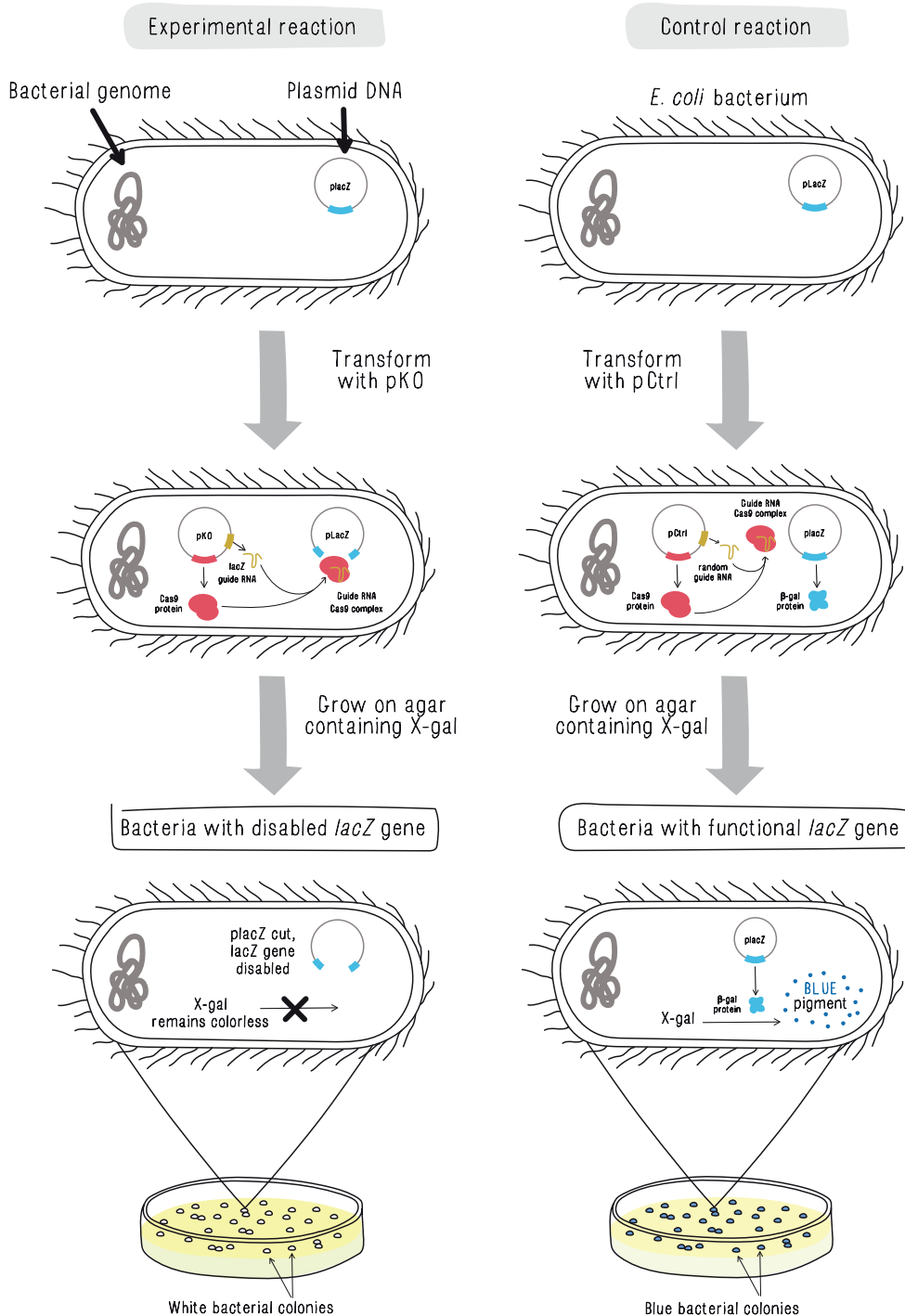


Figure 5. Experimental overview. Bacteria used in the experiment already contain the pLacZ plasmid. These cells are transformed with either pKO, which contains instructions for Cas9 protein and a guide RNA that targets *lacZ*, or pCtrl, which contains the instructions for Cas9 protein and a random guide RNA. The experimental transformation with pKO should lead to white colonies as the *lacZ* gene will be cut by Cas9 and disabled in the presence of a *lacZ* guide RNA (left). The control transformation with pCtrl should lead to blue colonies since the *lacZ* gene will remain functional in the presence of a random guide RNA (right).






Student lab protocol



Protective gloves and eyewear should be worn for the entirety of this experiment.

Part 1: Transformation

- You should receive two 1.5 ml tubes that contain bacteria. Label one tube “C” for control and the other tube “KO.”
- Place both tubes on ice.
- Add 25 μ l pCtrl DNA to tube “C.” There is no need to mix.
- Use a new pipette tip to add 25 μ l pKO DNA to tube “KO.” There is no need to mix.
- Securely close the lids of both tubes, then flick the tubes a few times to thoroughly mix the contents.
- Collect liquid in the bottom of the tubes by spinning briefly in a microcentrifuge or tapping the tubes on the tabletop.
- Incubate samples on ice for 30 minutes.
- Add 50 μ l recovery media to tube “C.”
- Use a new pipette tip to add 50 μ l recovery media to tube “KO.”
- Securely close the lids of both tubes, then flick the tubes a few times to thoroughly mix the contents.
- Collect liquid in the bottom of the tubes by spinning briefly in a microcentrifuge or tapping the tubes on the tabletop.

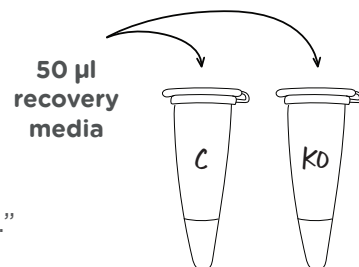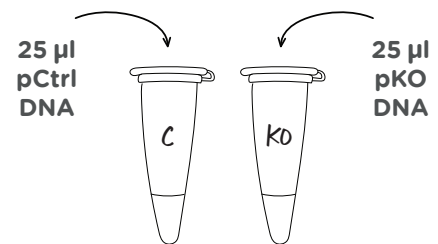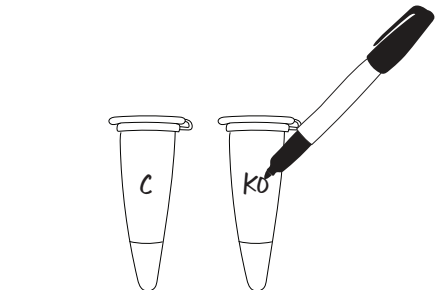

Step 12 of the protocol varies depending on whether you will plate your transformation reactions today or tomorrow. Ask your teacher which instructions to follow.

- If you will plate the transformation reactions today** (either during this class or at the end of the day): Incubate samples at 37 °C for at least 60 minutes. If you have access to a shaking incubator, you may shorten the incubation to 30 minutes.


OR

If you will plate the transformation reactions tomorrow: Leave the samples at room temperature on the lab bench overnight.



Part 2: Plate transformation reactions


13. You should receive two pre-warmed LB agar plates. The LB agar contains X-gal. Label the bottom of one agar plate “C” and the other agar plate “KO.” Also, label the bottom of each plate with your group name and the date.

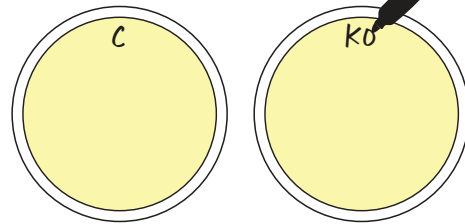

Tips for plating bacterial transformations


To open the spreader packaging:

- Peel back the clear part of the wrapper enough that you can grab the handle, but keep the head protected.

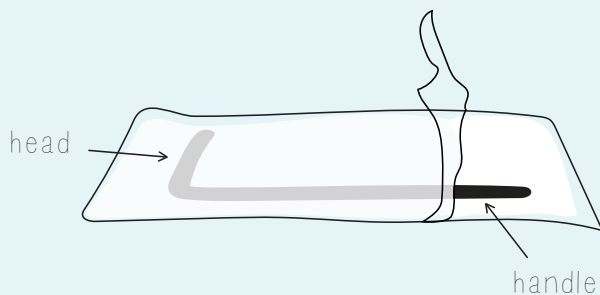

To use the spreader:


- Pipette transformation reaction onto the center of the plate.
- Lightly glide the head of the spreader across the agar—you do *not* need to apply pressure.
- Make sure to spread the transformation reaction over the entire surface of the agar.
- Spreaders are single use, discard them after one use.

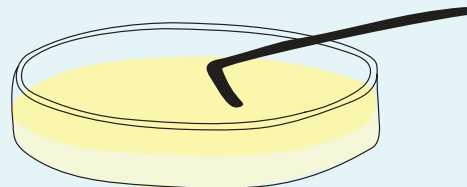

14. Plate control transformation reaction:

- Remove the lid from the plate labeled “C.”
- Use a micropipette to add the entire contents of tube “C” (~165 μ l) onto the agar plate labeled “C.”
- Spread mixture evenly across the surface of the agar using a sterile plastic spreader.
- Immediately place the lid back on the plate.
- Discard the spreader as instructed by your teacher.


~ 160 μ l control transformation reaction

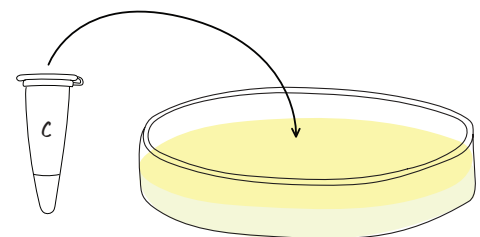

15. Plate KO transformation reaction:


- Remove the lid from the plate labeled “KO.”
- Use a micropipette to add the entire contents of tube “KO” (~165 μ l) onto the agar plate labeled “KO.”
- Spread mixture evenly across the surface of the agar using a sterile plastic spreader.
- Immediately place the lid back on the plate.
- Discard the spreader as instructed by your teacher.

~ 160 μ l KO transformation reaction

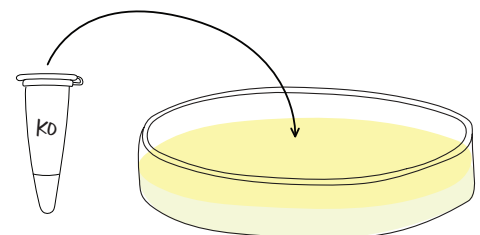



16. Let plates sit for at least 5 minutes to ensure that the liquid soaks into the agar.
17. Incubate plates upside down at 37 °C for at least 24 hours or until colonies are visible (this may take up to 48 hours).
 - Room temperature incubation is not recommended.
 - Once colonies are visible, plates can be stored in the refrigerator for several days.

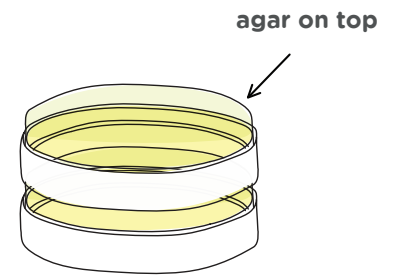

Part 3: Data collection

18. Record the approximate number of blue and white colonies on each plate in the table on page 31.
19. Dispose of plates as instructed by your teacher.



Study questions: pre-lab

Review

1. What is genome editing?

2. Why is the CRISPR/Cas genome editing system considered so revolutionary?

3. What are the two main components of the CRISPR/Cas genome editing system?

4. The CRISPR/Cas system has been compared to a homing missile programmed to find and damage a specific DNA sequence.

a. In this analogy, what part of the CRISPR/Cas system acts like the homing system, able to locate the target?

b. What part acts like the missile, damaging the target?



5. Why is it so important for scientists to target Cas9 with a high degree of specificity? What would happen if Cas9 cut somewhere other than at the intended target?

6. When Cas9 cuts DNA, the cell often introduces random mutations when DNA repair enzymes reattach the cut ends of the DNA back together. Explain why random mutations that inactivate a gene can still be useful to scientists.

7. When Cas9 cuts the DNA, this event is referred to as a “double strand break.” Why do you think this is the case?

8. We could use the CRISPR/Cas9 system to disrupt any gene, but this lab focuses on the *lacZ* gene because it is easy to know when the CRISPR/Cas9 system has worked effectively. Summarize how you can tell whether the *lacZ* gene is functional.



Critical thinking

9. Cas9 should not be able to cut DNA without a matching guide RNA. The cells that have Cas9 and random guide RNA are included in this experiment as a control.
- Why is it essential to include this control?

- Assume your results show white colonies after adding Cas9 and the random guide RNA. Brainstorm at least two possible explanations for this unexpected result.

Mathematical thinking

The product rule calculates the probability of a series of independent events by multiplying the probability of each event. If you flip a coin once, there is a 50% chance that it will land heads up. In fact, any time you flip a coin, the likelihood of getting a “heads” on that specific toss is 50%. If you want to know the probability of getting 3 “heads” in a row, you multiply $0.5 \times 0.5 \times 0.5$ to get 0.125 or 12.5%. The product rule can also be applied to calculate the probability of any DNA sequence. Because there are four possible DNA bases, the probability that any specific DNA sequence will occur is $(1/4)^n$, where n is the length of the sequence.

10. The part of the guide RNA that matches the target sequence and determines where Cas9 will cut the DNA is ~20 bases long. In the following questions you will explore how this allows scientists to target Cas9 to the desired regions of the genome with high specificity.
- What is the probability of any 20 base sequence occurring by chance? Show your work.

- The human genome is approximately 3.2 billion (3,200,000,000) bases long. How many times is any given 20 base sequence predicted to occur in the human genome? Show your work.



CER table

Fill in the table based on your results from the lab. Use the rubric on the next page to guide your answers.

Question:

Based on your results, on which plates did you successfully disrupt the *lacZ* gene?

Claim

Make a clear statement that answers the above question.

Evidence

Provide data from the lab that supports your claim.

Reasoning

Explain clearly why the data you presented supports your claim. Include the underlying scientific principles that link your evidence to your claim.



Score	4	3	2	1
CLAIM A statement that answers the original question/problem.	Makes a clear, accurate, and complete claim.	Makes an accurate and complete claim.	Makes an accurate but incomplete or vague claim.	Makes a claim that is inaccurate.
EVIDENCE Data from the experiment that supports the claim. Data must be relevant and sufficient to support the claim.	All of the evidence presented is highly relevant and clearly sufficient to support the claim.	Provides evidence that is relevant and sufficient to support the claim.	Provides relevant but insufficient evidence to support the claim. May include some non-relevant evidence.	Only provides evidence that does not support claim.
REASONING Explain why your evidence supports your claim. This must include scientific principles/knowledge that you have about the topic to show why the data counts as evidence.	Provides reasoning that clearly links the evidence to the claim. Relevant scientific principles are well integrated in the reasoning.	Provides reasoning that links the evidence to the claim. Relevant scientific principles are discussed.	Provides reasoning that links the evidence to the claim, but does not include relevant scientific principles or uses them incorrectly.	Provides reasoning that does not link the evidence to the claim. Does not include relevant scientific principles or uses them incorrectly.

We recommend that teachers use the following scale when assessing this assignment using the rubric. Teachers should feel free to adjust this scale to their expectations.

Rubric score	3	4	5	6	7	8	9	10	11	12
Equivalent Grade	55	60	65	70	75	80	85	90	95	100



Post-lab study questions

Interpreting results

1. Use the data table below to record your results

Plate	Approximate # of white colonies	Approximate # of blue colonies
Control rxn		
KO rxn		

2. Did you successfully disable the *lacZ* gene? Explain how you can tell.

Critical thinking

3. Were any of your results unexpected? If so, how were your results different from your expectations?



4. Use at least four of the following seven boxes to illustrate what occurred at the molecular level in this lab as if it were a comic strip. On the lines beside each box, describe what is happening in each drawing. Use and underline the following words or phrases:

- Cas9
- guide RNA targeting *lacZ*
- random guide RNA
- *lacZ* gene
- β -gal protein
- X-gal substrate
- blue colonies
- white colonies
- cut





Instructor guide



Expected results	P. 35
Troubleshooting	P. 36
Notes on lab design	P. 37
Additional student supports	P. 38
Extension activities	P. 38
Learning goals and skills developed	P. 39
Standards alignment	P. 39



Expected results

- Bacteria transformed with pCtrl (Cas9 and a random guide RNA) should form blue colonies
 - In the presence of a random guide RNA, Cas9 should not cut the *lacZ* gene.
 - The *lacZ* gene will remain intact and be used to express functional β -galactosidase protein.
 - In the presence of β -galactosidase protein, the X-gal substrate will be broken down, creating a blue product.
 - The blue product makes the *E. coli* colonies visibly blue.

- Bacteria transformed with pKO (Cas9 and a guide RNA targeting *lacZ*) should form white colonies
 - In the presence of a *lacZ* guide RNA, Cas9 should have cut the *lacZ* gene.
 - Cutting the *lacZ* gene disables it, so it can't be used to express β -galactosidase protein.
 - In the absence of β -galactosidase protein, the X-gal substrate will not be broken down, and X-gal will remain colorless.
 - *E. coli* colonies are naturally white.



pCtrl transformation



pKO transformation

Unexpected results

- Blue or partially blue colonies on pKO plates.
Because the *lacZ* gene is plasmid-based and the bacteria carry many copies of the pLacZ plasmid, it is possible that Cas9 will not cut and disable every copy of the *lacZ* gene. This partial knockout phenotype is more frequent when the cells are not incubated at 37°C, as well as if the cells are allowed to grow for more than 24 hours.

- White colonies on pCtrl plates
There is antibiotic selection for the pLacZ plasmid when the cells are rehydrated, but after the transformation, there is no selection for the pLacZ plasmid on the agar plates. Because the antibiotic selection on the agar plates is for the pCtrl or pKO plasmids, not the pLacZ plasmid, it is possible for the pLacZ plasmid to be lost from the cells once they are growing on the agar.



Troubleshooting

The transformation efficiencies for the plasmids used in this lab are lower compared to other plasmids. This is because the introduction of Cas9 can be somewhat toxic to the *E. coli* used in this experiment. Cas9 is a nuclease that cleaves DNA, and while it should only cut DNA that is complementary to the guide RNA used in the experiment, some off-target cutting can also occur. The *E. coli* strain used in this experiment lacks the ability to repair double-stranded DNA breaks, hence the cells will die if their genomic DNA is cut.

Because the transformation efficiencies are low, this experiment requires particularly precise adherence to the protocol. All of the following can adversely affect the transformation efficiency:

Initial growth period

It is essential for the bacteria to be in the mid-log phase of growth when used for transformation, and the protocol has been optimized to achieve this. Deviating from the recommended incubation of 16-24 hours at 37 °C can lead to transformation failure. It is also important that the cells be left undisturbed during this period as agitation accelerates growth.

Incubation on ice

Following the addition of DNA, it is essential that the transformation reactions stay ice cold. Using cubed ice instead of crushed ice or using pre-chilled metal blocks both lead to a notable decrease in transformation efficiency. Shortening the incubation on ice to less than 30 minutes also significantly decreases the transformation efficiency.

Recovery

Cells must be incubated for a short period under ideal growth conditions before plating. Allowing the cells to recover for at least 60 minutes at 37 °C in a standard incubator or at least 30 minutes at 37 °C in a shaking incubator increases transformation efficiency. Cells can also be left to recover overnight at room temperature before plating. Pre-warming the LB agar plates at 37 °C before plating noticeably improves transformation efficiency.



Notes on lab design

Focus on targeting: This lab emphasizes the use of guide RNAs to target Cas9 to specific DNA sequences. We have chosen to focus on targeting because it is central to the main advantage of the CRISPR/Cas system: its specificity.

Targeting a plasmid-based *lacZ* gene: Most lab strains of bacteria have been genetically modified so they do not repair DNA breaks. If the chromosomal DNA is cut, then these cells will die. Because of this, we chose to target a plasmid-based *lacZ* gene. The endogenous *lacZ* gene has been knocked out of the bacterial genome in many lab strains of bacteria, including the cells used in this lab. Thus, the *lacZ* gene is only found on the plasmid. In eukaryotic cells, after a double strand break the most common outcome is that non-homologous end joining (NHEJ) reattaches the cut ends of DNA. NHEJ typically introduces random mutations that inactivate the gene. Because lab strains of bacteria have been modified so they can't repair DNA breaks, after Cas9 introduces a double strand break in the pLacZ plasmid, the entire plasmid will be degraded.

Simplification of guide RNA structure: Throughout the lab and in the paper models, we represent the guide RNA as a single molecule. This is a simplification. In reality, the plasmids used for transformation contain the instructions for making two RNAs that work together. A CRISPR RNA (crRNA) serves to recognize the target sequence. The crRNA binds with another RNA called the transactivating CRISPR RNA (tracrRNA), which is responsible for interacting with Cas9.

IPTG induction of *lacZ*: In the pLacZ plasmid, the *lacZ* gene is inducible. The *lacZ* gene won't be transcribed unless the chemical IPTG is added. This inducible switch is based on the endogenous lac operon, and if you wish to discuss the regulation of gene expression with your students, you can tell them that IPTG has been added to the LB agar plates.

This lab was developed in partnership with Science Bridge, a former non-profit association based in Germany (<https://sciencebridge.net/>).



Additional student supports

E-worksheets: The student questions accompanying this lab are available for download [here](#) as editable text documents you can customize and upload to your LMS. E-worksheets can also be accessed from the Curriculum Downloads tab at <https://www.minipcr.com/product/knockout/>.

Bacterial transformation: A two-page primer for students who need a review of bacterial transformation and antibiotic selection. Available at <https://www.minipcr.com/tutorials/>.

CRISPR/Cas resource library: Available at <https://www.minipcr.com/CRISPR> to access the complete miniPCR bio™ CRISPR/Cas resource library, including:

- **Webinar: What is CRISPR?**
- **DNAdots article: CRISPR/Cas9**
Simple explanations of modern genetic technologies described in under two pages. Includes review, critical thinking, and discussion questions.

Extension activities

The following optional extension activities are provided for students to explore topics more deeply.

Paper models: We have created two paper activities for students to model and understand how CRISPR/Cas works. The paper models could be used before or after completing the lab. Download the paper models at <https://www.minipcr.com/crispr-paper-model/>.

- **CRISPR/Cas paper model: snail shell coiling:** To dive into a fun example of basic science research, students model the research mentioned in the background of the lab where scientists used CRISPR/Cas to determine the genetic control of the direction in which a snail's shell spirals (Abe and Kuroda, 2019). This activity is conceptually simpler than the next one.
- **CRISPR/Cas paper model: sickle cell gene therapy:** In terms of using CRISPR/Cas to treat a genetic disease, one frequently thinks of the more complex use of CRISPR/Cas to correct a mutation instead of using the system to knock out a gene. But scientists took the more basic knockout approach in the first patient to receive CRISPR/Cas gene editing to treat sickle cell disease. Students use the paper-based activity to model this cutting-edge therapy. The extension contains all the background information on sickle cell disease and hemoglobin for students to understand the mechanism underlying this gene therapy treatment for sickle cell disease. This activity is conceptually more advanced than the snail paper model as it introduces advanced aspects of gene regulation.

Knockout! PCR Genotyping Experiment (KT-1800-02): The optional PCR add-on provides a way for students to verify their phenotypic observations with molecular evidence.



Learning goals and skills developed

Student learning goals:

- Understand genome editing
- Describe why gene knockout is a powerful research tool
- Explain how CRISPR/Cas9 is used to target specific DNA sequences

Scientific inquiry skills:

- Students will identify dependent and independent variables and appropriate experimental controls
- Students will follow detailed experimental protocols
- Students will make a claim based in scientific evidence
- Students will use reasoning to justify a scientific claim

Molecular biology skills:

- Micropipetting
- Bacterial transformation
- Culturing bacteria
- Antibiotic selection
- Phenotypic screening

Standards alignment

The standards alignment document for this activity is available for download [here](#). This document can also be accessed from the Curriculum Downloads tab at <https://www.minipcr.com/product/knockout/>.

This activity is aligned to the following standards:

- Next Generation Science Standards: High School Life Science
- Advanced Placement Biology
- Texas Essential Knowledge and Skills: Biology
- Texas Essential Knowledge and Skills: Biotechnology
- Biotechnician Assistant Credentialing Exam
- Common Core ELA/Literacy Standards (9-10)

For additional information on alignment to state standards, please contact support@minipcr.com.